

Distribuido por:



Tel. +34 854 53 63 74
satcliente@lrdiagnostico.com
www.lrdiagnostico.com



CLIA-mate Plus

Chemiluminescence Immunoassay Analyzer

Aikang MedTech Co., Ltd



► CLIA-mate Plus

Principle	ALP-AMDDP chemiluminescence immunoassay
Sample type	Whole blood/Serum/Plasma
Throughput	30 t/hour
Test mode	Batch mode
Display	10.1-inch touch screen
Calibration frequency	100 days
Consumable	Disposable tips
Data storage	≥100,000
Communication port	USB/RS232/LAN/WIFI, available for LIS
Printer	Built-in thermal printer
Dimensions	575*325*605 mm
Weight	40 kg

► Easy Operation

Step 1



Prepare reagent cartridges
Put reagent cartridges into the sample tray

Step 2



Add sample
Add 80~100µl sample to the sample well of each reagent cartridge

Step 3




Insert reagent cartridges
Insert up to 6 reagent cartridges into the analyzer, add one pipette tip for each test

Step 4



Button start
Click "Start" on the touch screen, wait for the result

► All-in-one Reagent Cartridge

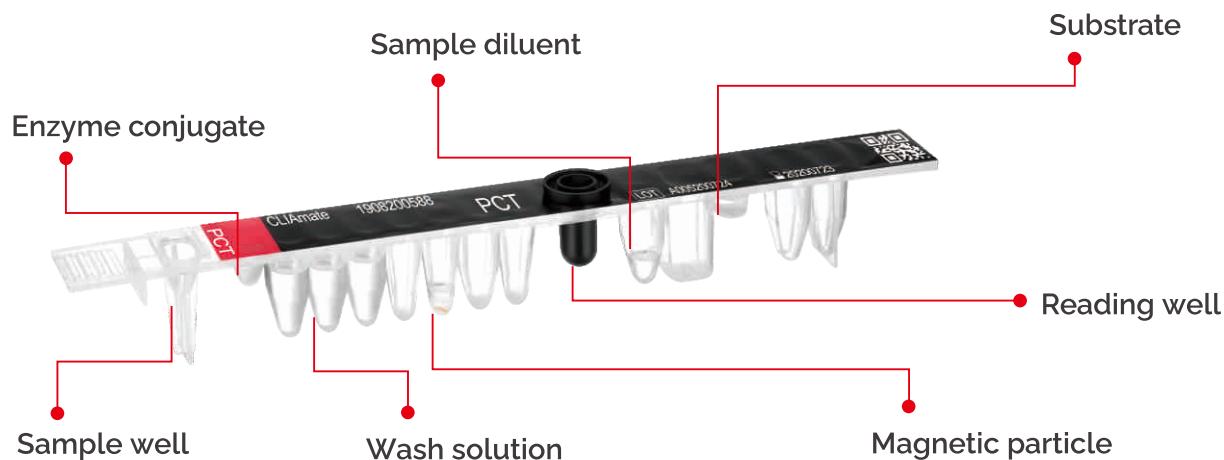


► Feature

- Real CLIA technology with high sensitivity&accuracy
- Instrument with compact design at affordable cost
- Single test reagent cartridge without any reagent waste
- Specification: 30 test/kit, 60 test/kit
- Disposable tip for each assay, reduce cross contamination
- Reagent information loading by external barcode scanner automatically
- 2 points calibration only, less calibration cost
- Android operation system, easy upgrading through USB port
- Easy operation and free maintenance
- Whole blood for cardiac and inflammation parameters, results within 17 minutes.



► Principle

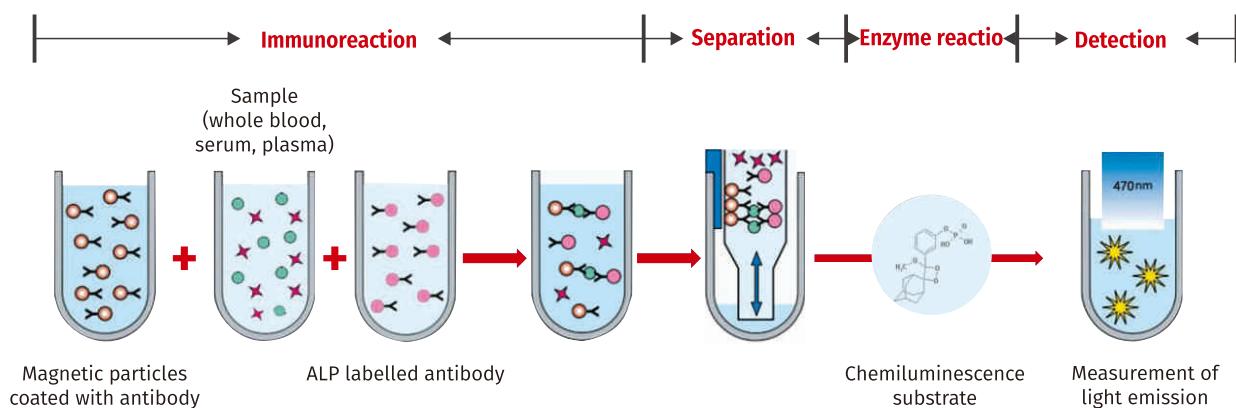

► Test Menu

► Thyroid

- TSH
- Total T3
- Total T4
- FT3
- Ft4

► Anemia

- VB12
- Folate*
- Ferritin

► Fertility

- AMH
- Progesterone
- β-HCG
- INH-B*
- PROL
- Testosterone
- FSH
- LH
- E2
- E3*

► TORCH

- TOXO IgG*
- TOXO IgM*
- RUB IgG*
- RUB IgM*

► Diabetes

- Insulin*
- C-Peptide*

► COVID-19

- IgG*
- IgM*
- Antigen*
- Neutralizing Antibody (nAb) *

► Infection Disease

- TB-IGRA
- HBs-Ag*
- Anti-HCV*
- HIV*
- TP*

► Tumor Marker

- CA125
- CEA
- tPSA
- fPSA
- AFP
- CA19-9*
- CA15-3*
- PG I
- PG II
- Gastrin-17

► Inflammation

- PCT
- CRP
- SAA
- IL-6
- HBP*

► Cardiac Markers

- cTnl
- Myo
- CK-MB
- NT-ProBNP
- D-Dimer
- cTnT
- H-FABP
- Hs-cTnl
- BNP
- Lp-PLA2*
- St2*
- MPO*
- Hcy*
- GDF-15*

► Brain Injury

- S100*
- Aβ1-42*
- AD7C-NTP*

► Autoimmune

- Anti-TPO*
- Anti-TG*

► Allergy

- IGE*

► Bone Metabolism

- 25-OH-VD
- BGP*
- PTH*

* Upcoming

Aikang Med Tech Co., Ltd.

- 📍 1A, Shijing Industrial Park, Pingkui Rd, Shijing St, Pingshan, 518118, Shenzhen, China
- 📞 +86 755 2790 1422
- ✉️ marketing@aikang-medical.com
- 🌐 www.aikang-medical.com

 AIKANG



All rights reserved. Features and specifications are subject to change without prior notice.
VN17.202204.02

Distribuido por:


Tel. +34 854 53 63 74
satcliente@lrdiagnostico.com
www.lrdiagnostico.com